


5 ml reaction tubes

The intermediate size for laboratories



The reaction tube for medium volumes:

- Easy handling of filling quantities between 2 and 5 ml
- Reliable SafeSeal closure
- PCR Performance Tested quality for more demanding uses
- Biosphere® Plus for maximum protection against contamination

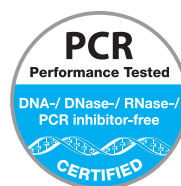

5 ml reaction tubes

Reliable handling of filling quantities between 2 and 5 ml

Routine applications in the laboratory such as dilution series, cell lysis or sample storing often require a filling quantity greater than 2 ml. Tubes with a larger volume are not often ideally suited for this due to their slim design or screw cap.

The new Sarstedt reaction tube combines the advantages of a SafeSeal reaction tube and a 5 ml screw cap tube in one. Easy to open with one hand thanks to the ergonomic design, the snap-on lids with a sealing lip and holding tab close securely even under thermal stress. The wide cap connector allows for easy and precise closing. Separate screw lids cannot be mixed up.



Advantages of the 5 ml reaction tube

- Robust polypropylene – can be centrifuged up to 21,000 x g, autoclavable
- Moulded graduation – easy to check volume
- Label field and frosted cap – safe labelling and quick identification
- Low height – low risk of contamination, the cone of the pipette does not have to be immersed in the tube
- Compact size – space-saving storage of samples, media or cultures

Order information

| Order no. | Name | Packaging |
|------------|--|------------------------|
| 72.701 | 5 ml reaction tube | 100 / bag 800 / box |
| 72.701.400 | 5 ml reaction tube, PCR Performance Tested | 100 / bag 800 / box |
| 95.1298 | 5 ml reaction tube centrifuge adapter | 8 / bag 96 / box |



Accessories

| Order no. | Name | Format | External dimensions (mm) | Packaging |
|-------------|---|--------|------------------------------|-----------|
| 93.850.001 | Polycarbonate rack for 18 x 5 ml reaction tubes | 6 x 3 | 137 x 70 x 40 (L x W x H) | 1 / box |
| 93.1431.001 | Polycarbonate rack for 48 x 5 ml reaction tubes | 12 x 4 | 257 x 90 x 40 (L x W x H) | 1 / box |

SARSTEDT AG & Co. KG
P.O. Box 12 20
D-51582 Nümbrecht
Phone: +49 2293 305 0
Fax: +49 2293 305 3992
export@sarstedt.com
www.sarstedt.com