


Low Profile PCR items

Optimized for 'Fast PCR'

NEW!

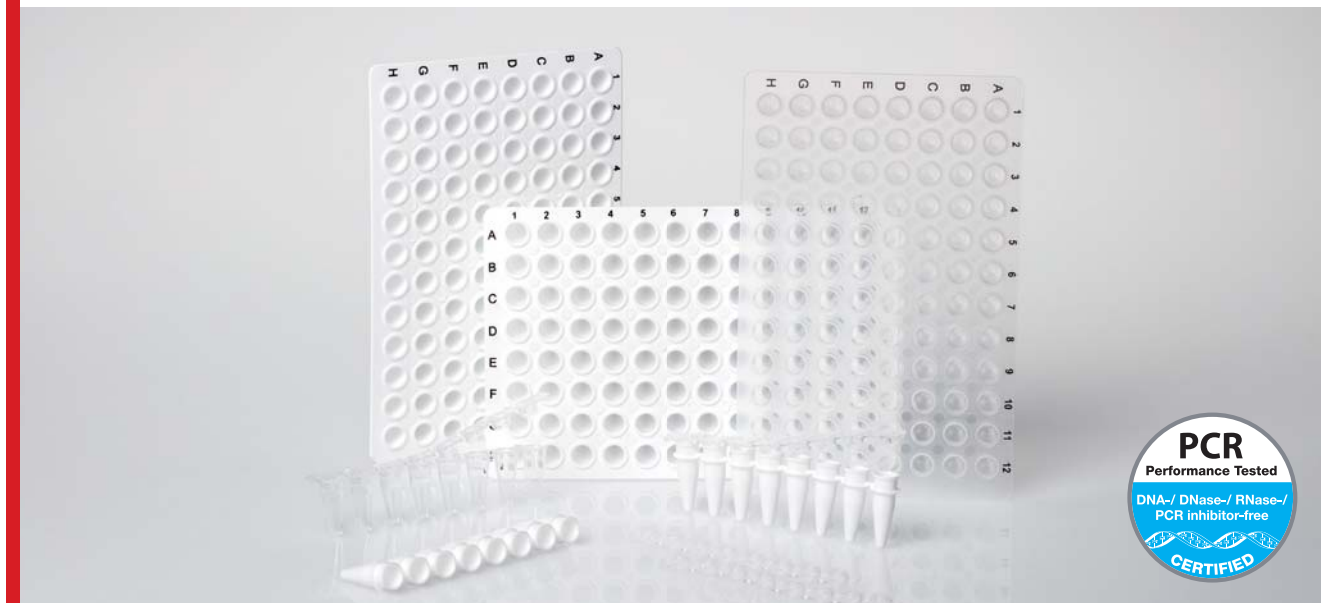

Sarstedt Low Profile PCR 8-strips and plates are specifically designed for small reaction volumes less than 20 μ l and are required for 'Fast PCR' thermocyclers (0.1 ml block format).

Product features:

- Due to the lower height, there is less air space available for the liquid PCR mixture to transition to the gas phase during heating, meaning that less condensation forms on the side walls and the reaction conditions are kept more constant.
- Extremely uniform and thin-walled design for a rapid temperature transition in the reaction space.
- For optimal fluorescence detection during qPCR there are white options available.
- Certified free of human DNA, bacterial DNA, DNase, RNase and PCR inhibitors specifically for the requirements for working with nucleic acids.

Ordering information

Order number	Description	Packaging
72.982.002	Multiply [®] - μ StripPro LP, transparent; lid strips included	125/1,250
72.982.092	Multiply [®] - μ StripPro LP, white; lid strips included	125/1,250
72.1977.202	Multiply [®] PCR Plate LP, transparent	20/100
72.1977.232	Multiply [®] PCR Plate LP, white	20/100

SARSTEDT, Inc.
1025 St. James Church Road
P.O. Box 468
Newton, NC 28658 - 0468
Toll free: +1 800 257 5101
Fax: +1 828 465 4003
customerservice.us@sarstedt.com
www.sarstedt.com

