


Urine Measuring System A06 W



User-friendly and precise urine monitoring system
providing optimum hygienic safety



Individual attachment

The perforated elastic straps and hook provide a safe and individual mounting of the system.

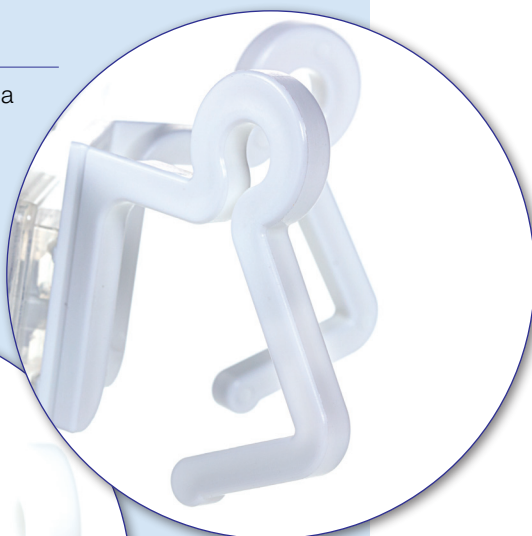

Hygienic & user-friendly exchange bag

The Bayonet-lock enables a quick and hygienic exchange of the bag and a secure connection to the measuring chamber.

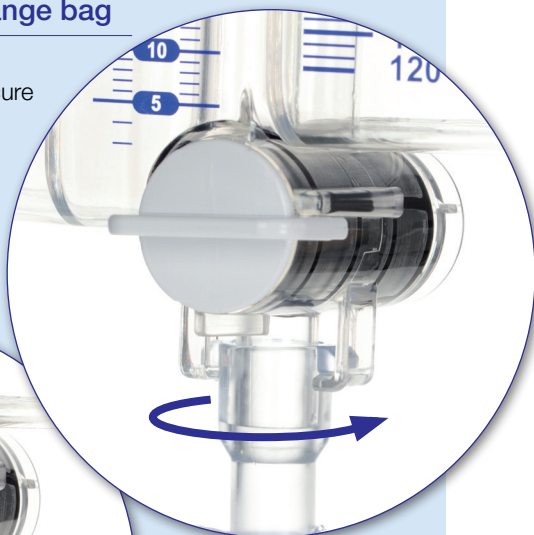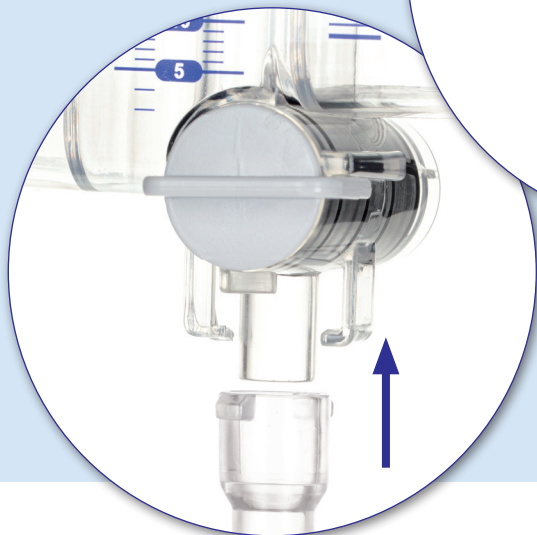



Needle-free sampling site

Needle-free, flat sampling site for secure disinfection and sampling.





Precise Reading

The precise and easily readable graduations facilitate accurate output readings.

The total measuring volume of 500 ml and the clear and concise graduations enable an accurate balancing of large urine volumes.






Volumes exceeding 500 ml are automatically diverted through the integral overflow chamber and drained into the collection bag.

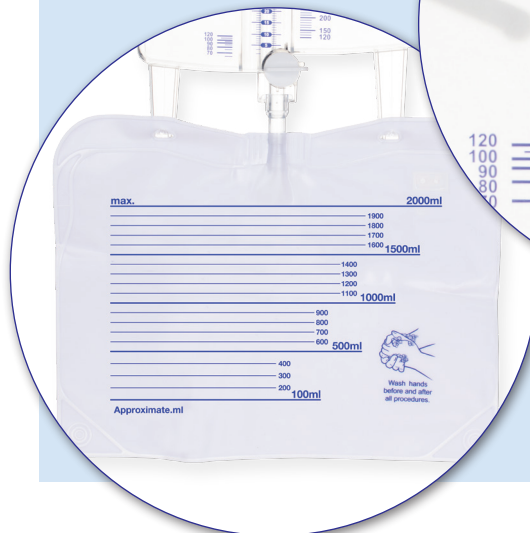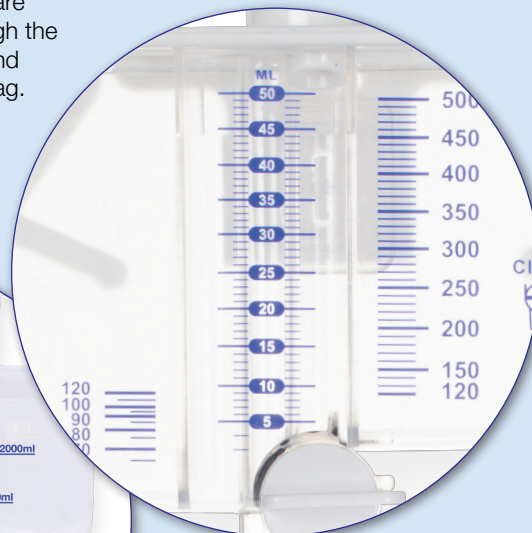

Special features

- Needle-free, flat sampling site
- Exchange bag is connected with a bayonet lock
- Universal mounting with elastic straps and additional hook
- Tube is kink resistant above the drip chamber
- Transparent drainage tap

Ordering Information

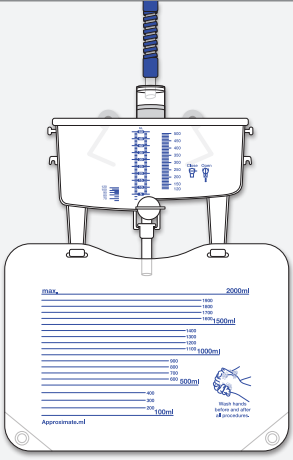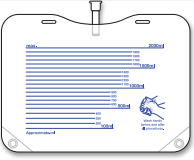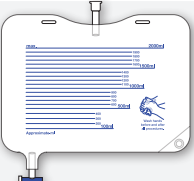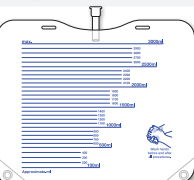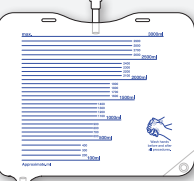	Order No.	Description	Packaging
	74.5681.200	Urine measuring system A06 W 500 ml with 2.0l exchange bag without drainage tap, with elastic straps and hook attachment	10
	74.5684.001	Exchange bag 2.0 l for A06 W without drainage tap	50
	74.5684.002	Exchange bag 2.0 l for A06W with drainage tap	50
	74.5685.001	Exchange bag 3.0 l for A06 W without drainage tap	50
	74.5685.002	Exchange bag 3.0 l for A06 W with drainage tap	50

Technical modifications reserved

This document may contain information on products that may not be available in particular countries

SARSTEDT AG & Co. KG
 P.O. Box 12 20
 D-51582 Nümbrecht
 Phone: +49 2293 305 0
 Fax: +49 2293 305 3992
 export@sarstedt.com
 www.sarstedt.com